


Dexcom G7 – How to share sensor data with the healthcare team

Using the Dexcom G7 phone app

Follow these steps to automatically share your Dexcom G7 data using the Dexcom G7 app:

1. Open the Dexcom G7 sensor app on your phone.
2. Tap the **Connections** tab.
3. Tap **Clarity Clinic** and tap **Continue**.
4. Enter the Clarity clinic code when prompted. The Best Centre code is: **bestdiabetes**
5. Confirm clinic details and tap **Done** when finished.

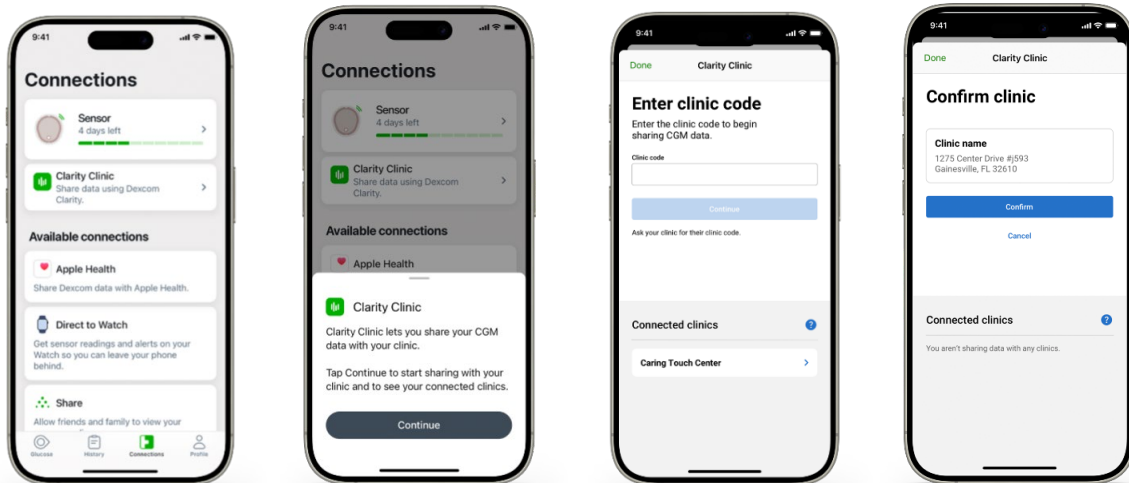

Important reminder:

If prompted to LOG IN to the app, always do so. You are not sharing data with your followers or the clinic if you are not logged in.

If you are having problems logging in, contact Dexcom Technical Support at 1-844-832-1810 for additional help.

Using the Dexcom G7 receiver

1. On a computer, go to <https://clarity.dexcom.eu/>. Select “Dexcom Clarity for Home Users”.
2. Log in to your account. If you do not have an account, click on “create an account” and follow the onscreen instructions
3. Click **upload** at the top of the screen.
4. If this is your first time using the Dexcom Clarity Uploader, click download to install it. Follow the onscreen instructions. When the installation is done, click on “refresh page”.
5. **Connect the receiver:** plug your G7 receiver into your computer using the USB charging cable that came with it.
6. Click **upload**, your receiver will automatically connect and send data.